



# APPLICATION NOTE

## MAESTRO FIRMWARE UPDATE PROCEDURE



ALWAYS KEEP  
YOUR MAESTRO  
UP TO DATE TO  
TAKE ADVANTAGE  
OF ALL THE NEW  
FUNCTIONALITIES

MONITORS

ENERGY DETECTORS

POWER DETECTORS

HIGH POWER DETECTORS

PHOTO DETECTORS

THz DETECTORS


OEM DETECTORS

SPECIAL PRODUCTS

BEAM DIAGNOSTICS

### FOLLOW THIS EASY STEP-BY-STEP UPDATE PROCEDURE

1. Download the update file from the website at:  
<http://gentec-eo.com/downloads/software-firmware>  
in the **Firmware** tab
2. Save it to a known location on your PC (like your desktop)
3. Unzip the file:
  - 3.1 Locate the ZIP file on your computer
  - 3.2 Right-click on it and choose "Extract All..."
  - 3.3 Follow the instructions on the screen
  - 3.4 Locate the extracted folder and open it
  - 3.5 You should now see the BIN file

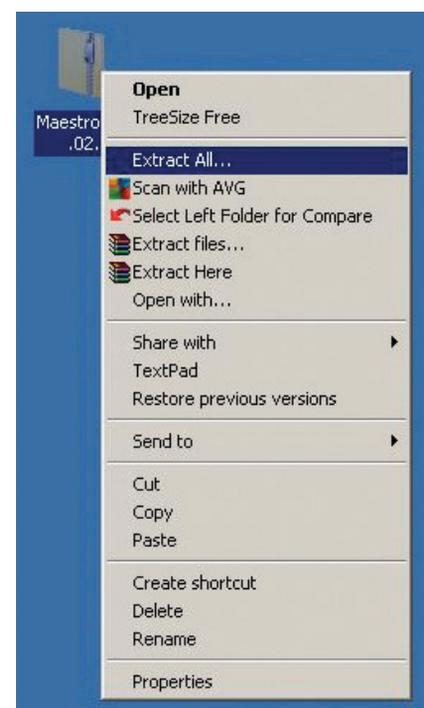



# APPLICATION NOTE

- Copy the BIN file on a USB key, directly in the main folder (not in a sub-folder)



- Plug the USB key in the MAESTRO

- Go to **HOME** 

- ABOUT** 

- Press the **Update** button 

- Wait for the MAESTRO to do the update sequence

- When the reboot message appears, turn off the MAESTRO

- Wait a few seconds

- Turn it on again and the update will be done

- If you want to check the update has been done, go in **HOME – About** again and verify that the firmware version corresponds to the initial update file